

# Quick Reference Guide

## Warning

### TO PREVENT INJURY:

READ THE OWNER'S OPERATING MANUAL BEFORE USING THE DENTAL AIR FORCE®.

PLACE THE DENTAL AIR FORCE® ON A SECURE, STRONG, FLAT SURFACE.

#### STEP 1



Screw the metal fitting on the coiled tubing into the metal fitting on the Dental Air Force®.

#### STEP 2



Twist the other end of the coiled tubing (with the plastic quick disconnect) onto the plastic fitting in the JetStream™ Tool Tip.

#### STEP 3




Fill the mixer ½ full with water.

#### STEP 4



Pour Dental Cleaning Agent into the Mixer  $\frac{1}{4}$  –  $\frac{1}{3}$  full. Fill any remaining space with water.

#### STEP 5

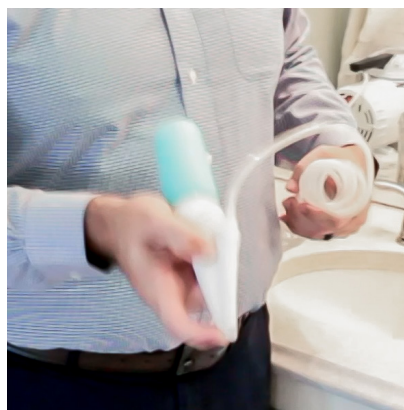

Assemble the JetStream™ tool by inserting the Tool Tip into the mixer and screwing it on securely. Plug the nozzle on the Tool Tip and shake thoroughly. Shake again during use.

#### STEP 6



Power on the Dental Air Force® and squeeze the mixer to control how much Dental Cleaning Agent is dispersed. Clean teeth for at least 2 minutes, twice a day.

